

CURIE Plus

**The Ultimate Biocompatible
UV Curing Oven for Dental 3D Printing**



Biocompatible Dental Workflows
Class I & Class II validated industry-leading resins automatically available for CURIE Plus



Chamber Capacity
Capable of curing 8 full arches simultaneously to maximize efficiency



Compatible with Nitrogen Sources
Built for nitrogen tanks and generators

Advanced Curing Technology

High-quality post-curing that ensures biocompatibility of Class II medical applications



65 LEDs with 360° light distribution



Mix-Matrix curing covering 365, 385, 405 nm UV wavelengths

A New Standard in Simplicity

CURIE Plus sets the standard for user experience while improving your productivity



Seamless integration with SOL for a 1-click workflow from print to cure



Cloud-based resin library and firmware updates

TECHNICAL SPECIFICATIONS



Dimensions

14.4 (W) x 10.2 (D) x 9.2 (H) in



Weight

26.9 lbs



Power Input

100-240VAC, 50-60Hz



Library Updates

Internet (Cloud)



UV Wavelength

365, 385, and 405 nm



Chamber Size

7.9 (W) x 5.9 (D) x 2.4 (H) in
Fits up to 8 full arches




Technology

UV LED Array

CURIE Plus - Item Number

RS110010162

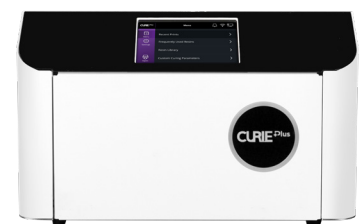

PACKAGE CONTENTS

- 1 x Glass Panel Box
- 1 x Power Cable
- 1 x Ethernet Cable
- 1 x Wi-Fi Dongle
- 1 x N2 Cover Kit



Scan the QR Code to
Schedule a Virtual Demo Today!



Ackuretta