### Stabilimento

CEFLA NORTH AMERICA Inc, 6125 Harris Technology Blvd., 28269 Charlotte NC Tel. +1 (704)598-0020 Sales Support +1(800)416-3078

#### **BU Medical Equipment** Sede Legale Ed Amministrativa Headquarters

Tech Service +1(866)466-5228

CEFLA s.c. - Via Selice Provinciale, 23/a - 40026 Imola (B0) - Italy Tel. +39 0542 653111 Fax +39 0542 653344


#### Stabilimento

CEFLA s.c. - Via Bicocca, 14/C - 40026 Imola (BO) - Italy Tel. +39 0542 653441 Fax +39 0542 653601

#### CASTELLINI.COM

# PUMA ELI AMBIDEXTROUS

PERFECT POSITIONING FOR ANY STYLE

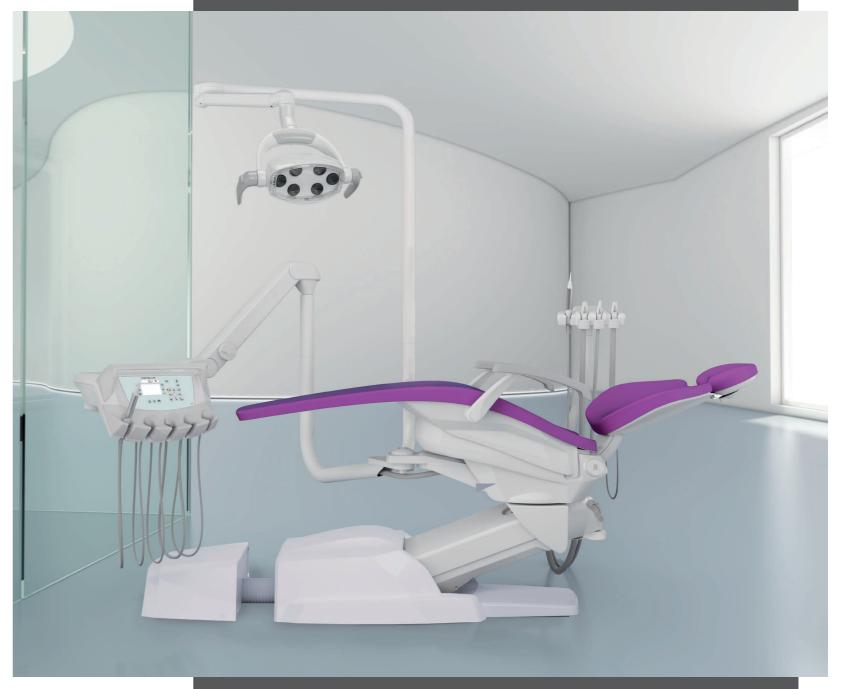





## BALANCE AND COMFORT

Compact, light and practical - also in the switchover from right-handed to left-handed user configuration - the Puma ELI Ambidextrous dental unit retains an essential streamlined style and ensures performance always meets the dentist's every clinical requirement.

Free to change.







#### ADAPT TO YOUR LAYOUT

Free to position the versatile Junction Box anywhere in the operatory, dental surgeons can adapt easily to the layout and connect the unit to all chair utilities as required. Such freedom represents a simple solution today and avoids limitations if the layout changes in the future.

The Junction Box can be placed at up to nearly 60 inches from the base of the dental unit.



#### PIVOT

The patient chair can rotate 30° to the left or right, making access easier for patients. It also helps the dental team to adjust the position according to the required treatment.



#### S.S.S.

The device that feeds water to the sprays as an independent alternative to mains water is supplied as standard. Works by way of a tank filled with distilled water: this prevents limescale build-up. Extremely useful where mains water is hard.



#### STANDARD HEADREST

The easy-to-adjust headrest can be adapted to patients of different heights.



#### KEYPAD

Located on the side of the seat, this keypad lets users control patient chair movement and allows set positions to be saved and recalled.

PUMA ELI AMBIDEXTROUS PERFECT POSITIONING FOR ANY STYLE



#### TOUCHSCREEN CONTROL PANEL

Puma ELI Ambidextrous is equipped with a new instrument control panel with touchscreen and integrated X-ray viewer. An alphanumeric display allows for smooth conservative dentistry function control, providing essential data and allowing operating parameters to be set directly. The "clean" button temporarily disables the touchscreen to allow sanitization.



#### ASSISTANT'S MODULE

The assistant can manage a wide range of functions via the onmodule control panel. This lets dentists focus their attention on actual treatment. In addition to patient chair movement controls, there is also controls for the lighting,



#### C8 ST00

The C8 stool is designed to provide maximum mobility, quick access to all instruments and optimal arm support.

#### C7 STOOL

The C7 is an optimized stool that lets dentists perform treatment perfectly while enjoying correct body support. Adjustment mechanisms ensure both seat and back are always regulated to suit the operator's build, while the compact design allows full access to all parts of the dental unit. The lining material and the manufacturing technique make it easy to sanitize the surfaces.



#### POWER PEDAL FOOT CONTROL

The power pedal foot control is characterized by two additional elements to recall two more saved working positions.



#### VENUS E LED OPERATING LIGHT

Supplied as standard with 3-axis movement, front panel, removable handles and six LED sources to minimize shadows. Light intensity can also be set, via the sensor, within the 7,000 - 40,000 Lux range. The No Touch sensor lets the user switch the lamp on/off without touching any button.



#### SYRINGE

Featuring a streamlined and easy-to-handle design, the syringe is available both in the three-way and the optional six-way version. The handle and tip are autoclavable.



#### MEMORY FOAM

Available as an optional, Memory Foam padding adapts to patients of any build, ensuring proper anatomical support and maximum comfort at all times.



102 Atlantic blue 113 Pacific blue

106 Mediterranean blue136 Indian blue

101 Caribbean green 123 Polynesian green

123 Polynesian green103 Nevada yellow

115 Scottish salmon 132 Blueberry violet

132 Blueberry violet134 Japanese wisteria135 Venetian red

130 Graphite black121 Anthracite grey137 Satin silver

MEMORY FOAM PADDING

198 Atlantic blue 183 Pacific blue 196 Mediterranean blue

186 Indian blue197 Caribbean green193 Polynesian green

182 Nevada yellow

195 Scottish salmon

192 Blueberry violet 184 Japanese wisteria

194 Venetian red

180 Graphite black 199 Anthracite grey

187 Satin silver



PUMA ELI AMBIDEXTROUS PERFECT POSITIONING FOR ANY STYLE